



PMC-150

Cooled High Speed PMT Module for TCSPC

Fast TCSPC Instrument Response: < 150 ps FWHM

Internal Cooler: Low Dark Count Rate

Internal GHz Preamplifier: High Output Amplitude

Internal High Voltage Power Supply


Direct Interfacing to all bh Photon Counting Devices

Power supply, Gain Control and Overload Shutdown via bh DCC-100 module

Standard C Mount Adapter



The PMC-150 is a cooled PMT module for TCSPC applications. It contains a fast miniature PMT along with a Peltier cooler, a high voltage generator, a GHz pulse amplifier and a current sensing circuit. Due to the high gain and bandwidth of the device a single photon yields an output pulse with an amplitude in the range of 100 to 200 mV and a pulse width of 1.5 ns. Due to the high gain and the efficient shielding noise pickup is minimised. Therefore the PMC-150 yields high time resolution and high counting efficiency. The TCSPC instrument response function (IRF) has a width of less than 150 ps FWHM. Overload conditions are detected by sensing the PMT output current. Overload is indicated by an LED, an acoustic signal, and a logical overload signal. The PMC-150 is operated by the bh DCC-100 detector controller card. The DCC delivers the operating voltage for the PMT, the current for the Peltier cooler, controls the detector gain, and shuts down the PMT in case of overload. Compared to its predecessor, the PMC-100, the PMC-150 has a shorter IRF width and a better IRF uniformity over the active area.




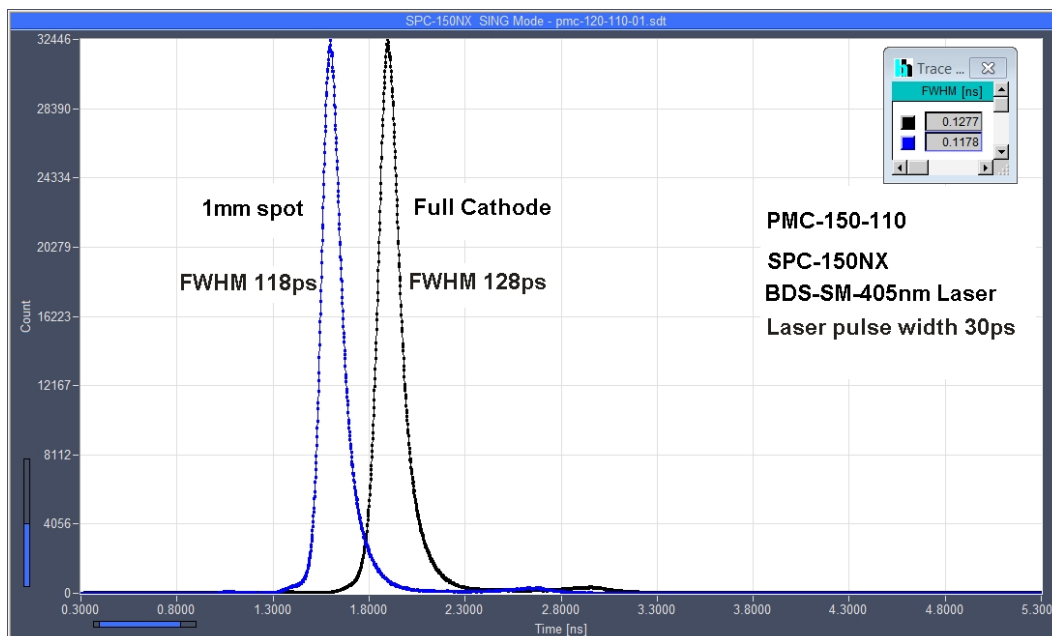


Becker & Hickl GmbH
Nunsdorfer Ring 7-9
12277 Berlin
Tel. +49 / 30 / 212 800 20
FAX +49 / 30 / 212 800 213
<http://www.becker-hickl.com>
email: info@becker-hickl.com



PMC-150

	PMC-150-04	PMC-150-110	PMC-150-113	PMC-150-210	PMC-150-20
Cathode type	Multialkali, UV	Super Bialkali	Super Bialkali, UV	Ultra Bialkali	Extended Red
Wavelength Range (nm)	185 to 870	230-700	185 to 700	230 to 700	300 to 900
Dark Counts (Icool = 1A, Tamb = 22°C, typ. value)	20	5	5	10	500
Cathode Diameter			8 mm		
IRF width, 1mm spot / full cathode area, typ., FWHM			120 ps / 130 ps		
Transit time variation with x-y position, see diagrams			50 ps		
Single Electron Response Width			1.5 ns, FWHM, typ. value		
Single Electron Response Amplitude			100 to 200 mV, Vgain = 0.9 to 1V		
Output Polarity			negative		
Count Rate (Continuous)			> 5 MHz		
Count Rate (Peak, < 1 us)			> 100 MHz		
Overload Indicator			LED and acoustic signal		
Overload Signal			TTL / CMOS, active low		
Detector Signal Output Connector			SMA		
Output Impedance			50 Ω		
Power Supply (from DCC-100 Card)			+ 12 V, -12V (fan only), Cooler current 0.5 to 1A		
Dimensions (width x height x depth)			76 mm x 111 mm x 56 mm		
Optical Adapter			C-Mount female		
Fibre Coupling			SMA 905 or FC, on request		

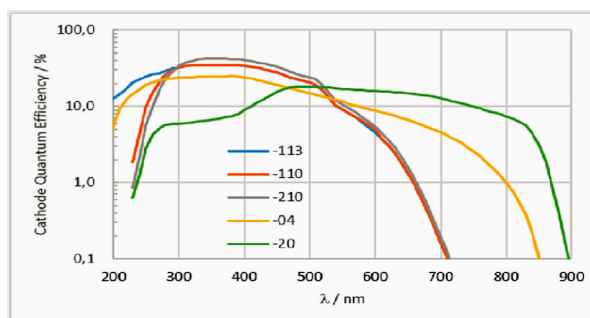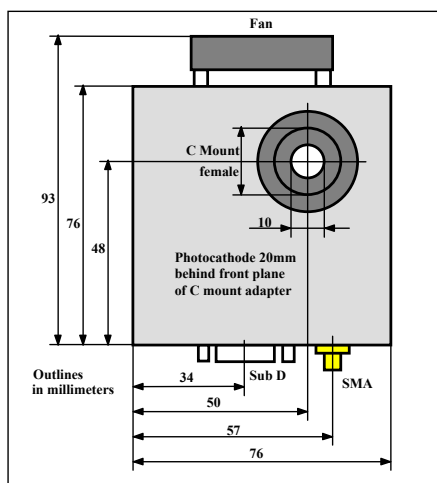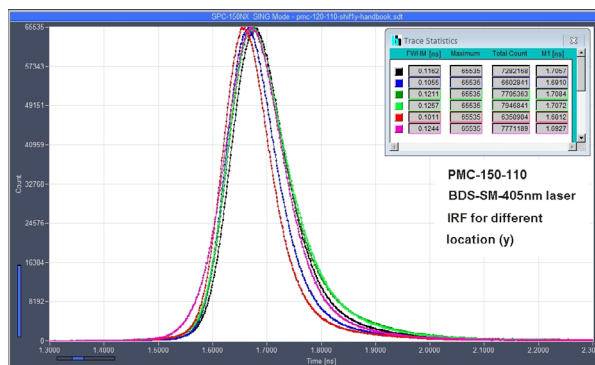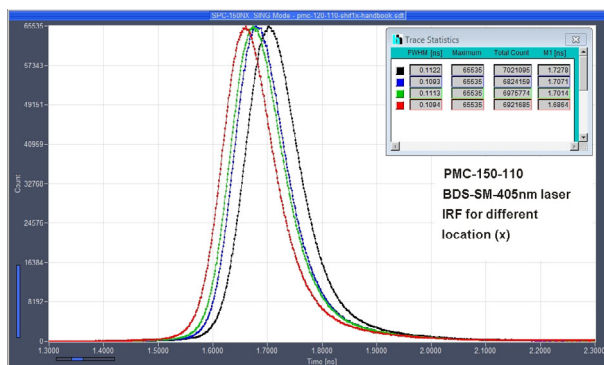

Pin Assignment of 15 pin sub-d-hd connector

1	not used	9	Peltier -
2	Peltier +	10	+12V
3	Peltier +	11	-12 (Fan)
4	Peltier +	12	not used
5	GND	13	Gain Control, 0 to +0.9V
6	not used	14	/OVLD
7	Peltier -	15	GND
8	Peltier -		

International Sales Representatives



US:
Boston Electronics Corp
tcspe@boselec.com
www.boselec.com



UK:
Photonic Solutions PLC
sales@psplc.com
www.pspcl.com



TII Group Company



UNISOKU 超高真空・極低温走査型プローブ顕微鏡
高速分光測定装置、クライオスタット



LOTIS TII Nd:YAGレーザー、Ti:Sレーザー
OPOレーザー

東京インスツルメンツ TOKYO INSTRUMENTS 本社: 〒134-0088 東京都江戸川区西葛西6-18-14T, ビル Tel. 03-3686-4711
営業所: 〒532-0003 大阪府大阪市淀川区宮原4-1-46 新大阪北ビル Tel. 06-6393-7411
グローバルにネットワークを広げ、最先端の科学をお客様に提供 URL: https://www.tokyoinst.co.jp Mail: sales@tokyoinst.co.jp

●本カタログに記載されている内容は、改良のため予告無く変更する場合があります。(製品の仕様、性能、価格などはカタログ発行当時のものです)
●本カタログに記載されている内容の一部または全部を無断で転載することは禁止されています。
●本カタログに記載されているメーカー名、製品名などは各社の商標または登録商標です。