


Time Tagger Series

Streaming time-to-digital converters

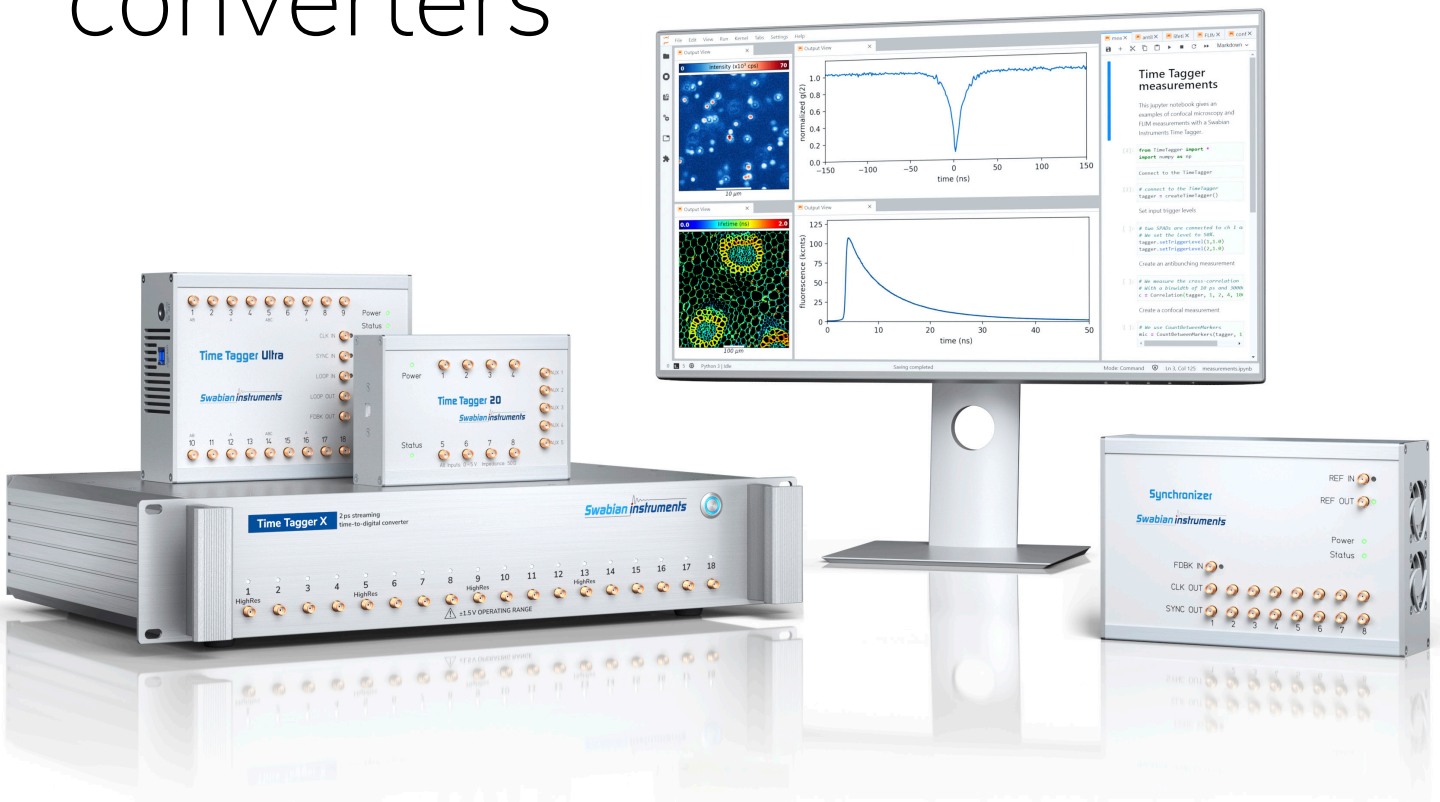

TII 東京インストルメンツ
TOKYO INSTRUMENTS
 グローバルにネットワークを広げ、最先端の科学をお客様に提供

本 社：〒134-0088 東京都江戸川区西葛西6-18-14 T.ビル Tel. 03-3686-4711
 営業所：〒532-0003 大阪府大阪市淀川区宮原4-1-46 新大阪北ビル Tel. 06-6393-7411
 URL: <https://www.tokyoinst.co.jp> Mail: sales@tokyoinst.co.jp

TII Group Company

UNISOKU 超高真空・極低温走査型プローブ顕微鏡
 高速分光測定装置、クライオスタット

LOTIS TII Nd:YAGレーザー、Ti:Sレーザー
 OPOレーザー


- 本カタログに記載されている内容は、改良のため予告無く変更する場合があります。(製品の仕様、性能、価格などはカタログ発行当時のものです)
- 本カタログに記載されている内容の一部または全部を無断で転載することは禁止されています。
- 本カタログに記載されているメーカー名、製品名などは各社の商標または登録商標です。

TIME TAGGER SERIES



Streaming time-to-digital converters


Swabian Instruments' Time Taggers are streaming time-to-digital converters with a unique data processing architecture that makes them the preferred choice for Time-Correlated Single-Photon Counting (TCSPC), time-interval counting, coincidence counting, and digital protocol analysis.

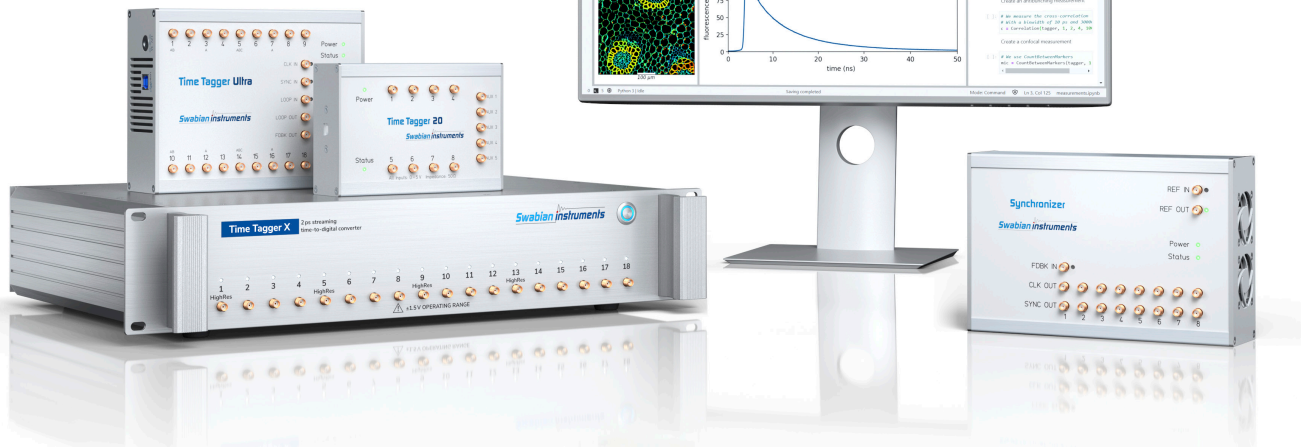

The new standard in Time-Correlated Single Photon Counting.

The Time Tagger Series provides endless capabilities for single-photon counting and you unleash them with no efforts. Whether you use the Time Taggers' powerful Software or derive from extensive code examples in Python, Matlab, LabVIEW, or C#/C++ - you get your experiments up and running within minutes.


- » antibunching
- » coincidence counting
- » laser scanning microscopy
- » multidimensional histograms
- » supports SNSPDs, APDs, SPADs, PMTs, SiPMs

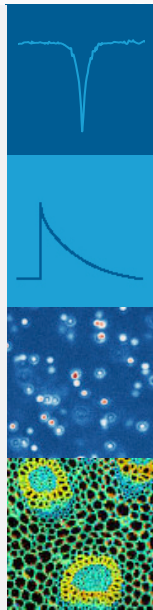

Timing & Synchronization made easy.

Are you looking for a versatile frequency counter or continuous time interval analyzer (CTIA)? Do you perform Allan Deviation, Modified Allan Deviation, or Hadamard Deviation measurements? Or do you test synchronization of 1 PPS signals? The Time Tagger Series lets you do this easier than ever.


- » ADEV 4E-13 at 1 s, TDEV 1E-13 at 1 ms
- » realtime ADEV, MDEV, HDEV, phase and frequency error measurements
- » up to 160 inputs simultaneously
- » input frequencies from 1 Hz to 700 MHz
- » easy 1 PPS testing and logging
- » log and post process your raw data

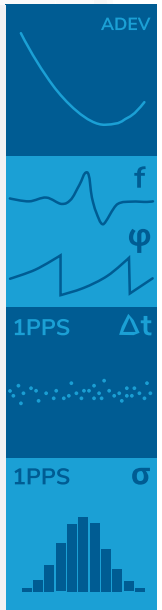

Take Dynamic Light Scattering to the next level.

Unsurpassed timing resolution paired with unlimited logarithmic histograms on a multitude of detectors. These are the Time Tagger's capabilities for dynamic light scattering.

- » measure simultaneously at up to 18 different scattering angles
- » remove intensity spikes on-the-fly and in post processing with advanced data filtering tools
- » calculate particle size distributions with Cumulants and CONTIN analysis
- » resolve down to picoseconds

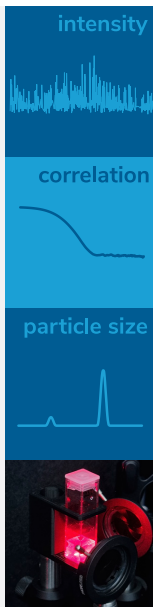

Endless processing capabilities.

Learn more at

swabianinstruments.com/static/documentation/TimeTagger/

Time Tagger 20



Time Tagger Ultra



Time Tagger X



Synchronizer



Compact. USB powered.

34 ps
RMS jitter

9 M tags/s
over USB 2.0

8
input channels

Flexible. Powerful.

Value

Performance

42 ps
RMS jitter

90 M tags/s
over USB 3.0

min **4** and up to **18**
input channels

Over-the-air upgradable

+ High Resolution (3 ps)

+ more input channels

+ Value to Performance

X-treme resolution. X-treme interfaces.

2 ps
RMS jitter

90 M tags/s
over USB 3.0

min **4** and up to **20**
input channels

Over-the-air upgradable

+ High Resolution (1.5 ps)

+ more input channels

Large scale systems. Combined into one.

0.8 ppb
clock stability over 24 h

combine up to **8**
Time Tagger Ultra and X

up to **160**
input channels total

+ rack mount kit



Implement your ideas within minutes

What makes our Time Taggers unique is their powerful software engine that offers effortless data processing capabilities. You will implement your measurement ideas within minutes - promised.



High timing resolution

The low instrument jitter down to 1.5 ps RMS (3.5 ps FWHM) combined with exceptionally short dead time down to 1.5 ns makes sure that you achieve the highest timing resolution for your application.



Versatile onboard event filter

A unique onboard event filter enables you to increase the effective time tag rate way beyond the streaming bandwidth by filtering out unneeded time tags right on the hardware.



Unlimited network capabilities

You can use the Time Tagger's software engine to cast your Time Tag Streams into a network. On client computers you can fire up virtual Time Taggers and work with them as if they were physically attached to your computer, providing the full measurement and data processing capabilities of a hardware Time Tagger.



Maximum flexibility

Run measurements independently using any combination of your input channels. You can take data simultaneously from independent physical setups and run multiple measurements on the same channels at the same time.



High data transfer rate

The high data transfer rate of 90 M tags/s to your computer over USB 3.0 enables you to process huge amounts of events on-the-fly. Leverage our SFP+ and QSFP+ interfaces for low latency experiments.



Native software libraries

Run your measurements in your preferred programming language with our included native software libraries, covering Python, MATLAB, LabVIEW, C#, C++, and even Mathematica.



Low latency FPGA output

For fast-feedback, low latency (< 10 μ s) experiments, transfer time tags to your own FPGA and enable endless processing capabilities via our SFP+ (300 M tags/s) and QSFP+ data interfaces (1200 M tags/s). Check our FPGA reference designs to get started immediately.

Timing precision	Time Tagger 20	Time Tagger Ultra		Time Tagger X
RMS jitter (typical)	34 ps	42 ps (Value)	8 ps (Performance)	2 ps
RMS jitter (typical, HighRes)	-	3 / 4 / 6 ps (2 / 4 / 8 HighRes channels)		1.5 ps (5 HighRes channels)
FWHM jitter (typical)	80 ps	100 ps (Value)	19 ps (Performance)	4.7 ps
FWHM jitter (typical, HighRes)	-	7 / 10 / 14 ps (2 / 4 / 8 HighRes channels)		3.5 ps (5 HighRes channels)
digital resolution	1 ps	1 ps		1 ps

Processing capabilities

input channels	8	4 to 18		4 to 20
dead time	6 ns	2.1 ns		1.5 ns
data transfer rate (to PC)	9 M tags/s	90 M tags/s		90 M tags/s
data transfer rate (SFP+, QSFP+)	-	-		300, 1200 M tags/s
burst memory	8 M tags	512 M tags		512 M tags
maximum input frequency	167 MHz	475 MHz		700 MHz

Input signals

input impedance	50 Ω	50 Ω	50 Ω / 1 MΩ
recommended input signal range	0 to 3 V	-3 to 3 V	-1.5 to 1.5 V
input signal range	-0.3 to 5 V	-5 to 5 V	-3 to 3 V
trigger level range	0 to 2.5 V	-2.5 to 2.5 V	-1 to 1 V
minimum pulse width	1000 ps	500 ps	350 ps
minimum pulse height	100 mV	100 mV	100 mV

External clock input

frequency	-	10 MHz or 500 MHz	10 MHz or 500 MHz
coupling	-	AC, 50 Ω	AC, 50 Ω
amplitude	-	1 to 3 Vpp	0.5 to 4 Vpp







General parameters

data interface	USB 2.0	USB 3.0	USB 3.0, SFP+, QSFP+
----------------	---------	---------	----------------------

Typical performance

Instrument response

1 MHz square wave, 1 Vpp, 1 ns rise, applied to two input channels, trigger 50%. The standard deviation σ of the distribution measures the jitter of two input channels. The RMS jitter of each individual channel is $\sigma/\sqrt{2}$. The FWHM jitter of each channel is $2.35 \sigma/\sqrt{2}$.

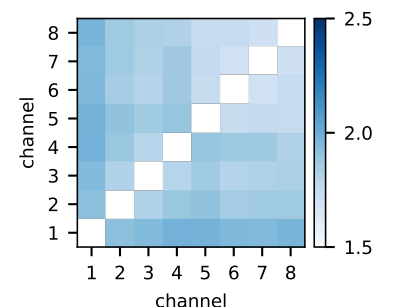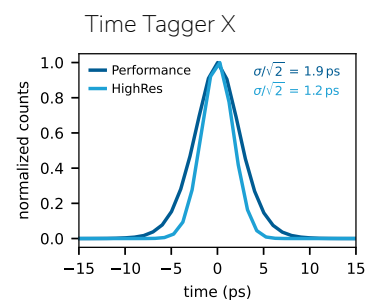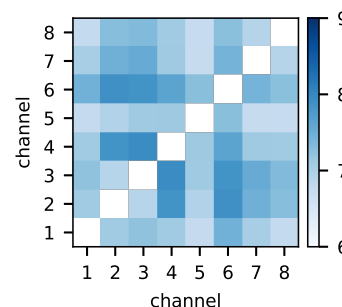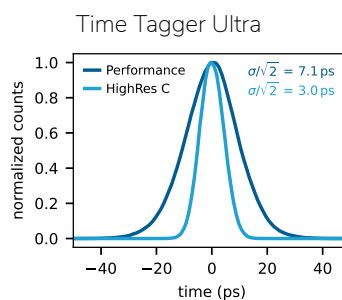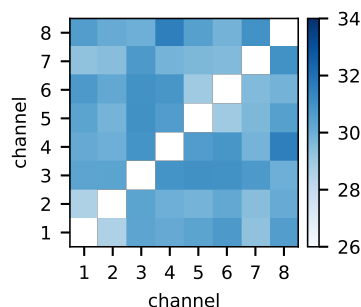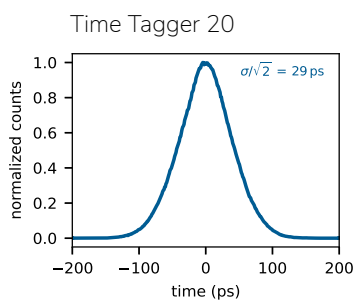

RMS Jitter

The plots show the RMS jitter obtained from instrument response measurements with all pairs of the first 8 input channels.